


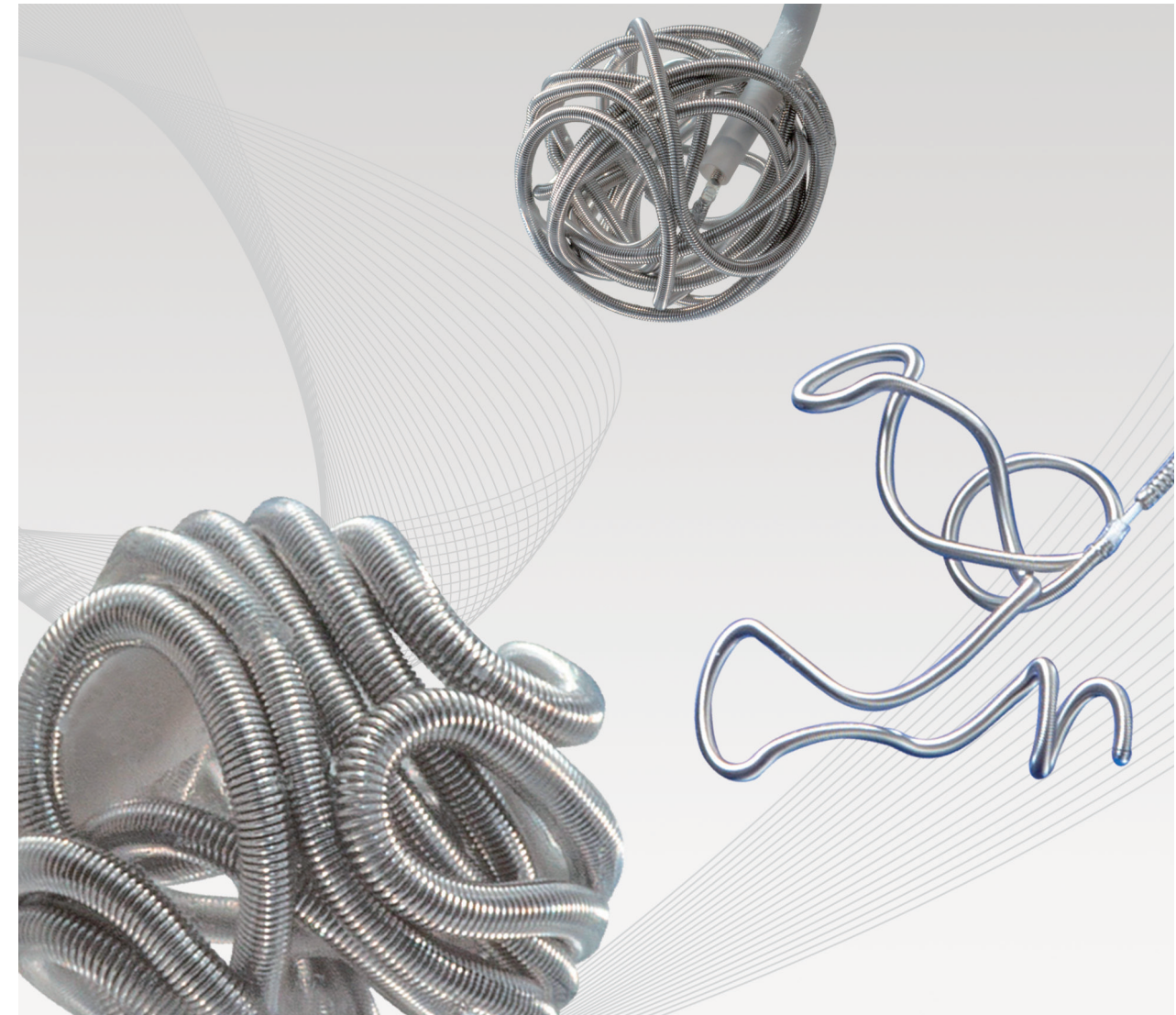


Product Specifications

✓ Pusher Length 187cm




Type	Catalogue No.	Primary Coil Diameter (inch)	Secondary Coil Diameter* (mm)	Primary Coil Length (cm)
 Complex	ExtraSoft (ES)	0.010	3.5	392-0304: 4
	392-0306: 6			
	392-03H8: 8			
	392-0408: 8			
	392-0410: 10			
	392-0510: 10			
	Soft (S)	0.012	5	392-0515: 15
	392-0620: 20			
	392-0730: 30			
	392-0840: 40			
	392-0940: 40			
	392-1040: 40			
MediumSoft (MS)	0.014	10	392-1250: 50	
Standard (STD)			16	392-1450: 50
				392-1650: 50
 Complex	SilkySoft (SS)	0.010	2.5	390-0101: 1
	390-0102: 2			
	390-0103: 3			
	390-01H2: 2			
	390-01H3: 3			
	390-0202: 2			
	390-0203: 3			
	390-0204: 4			
	390-0206: 6			
	390-0208: 8			
	390-02H4: 4			
	390-02H6: 6			
	390-0304: 4			
	390-0306: 6			
	390-0308: 8			
	 Complex Infini			SilkySoft (SS)
393-020306: 6				
393-020310: 10				
393-030515: 15				
393-030520: 20				

* The secondary coil diameter means referenced diameter of a space to be embolized.



i-ED COIL™
Detachable Coil

Type	Catalogue No.	Description
EDG v4	347-404	ELECTRO DETACH GENERATOR v4

 Do not re-use
  Consult instructions for use
  Sterilized using ethylene oxide

Wholesaler Vasnova
 82, Bugok-ro, Pogok-eup, Cheoin-gu,
 Yongin-si, Gyeonggi-do, Korea (PO) 17023
 E-mail. info@vasnova.co.kr
 Tel. +82.02.6959.5858 | Fax. +82.02.2038.3011

Manufacturer KANEKA CORPORATION
 3-18, 2-Chome, Nakanoshima, Kita-ku, Osaka-city,
 OSAKA, 530-8288 JAPAN
 Tel. +81.6.6226.5256
<https://www.kaneka-med.jp/english/>



Complex ExtraSoft / Soft / MediumSoft / Standard




Framing that is stable & easily fits into aneurysms

- Designed to sit in open space
- The distal coil shape spreads along the aneurysm wall & the proximal coil shape provides support from the inside.
- Coil designs optimized according to size
- Combination of Primary Coil Diameter & length for easy delivery & stability.



Complex SilkySoft

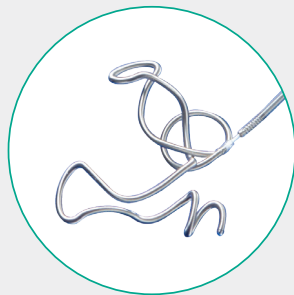

Pursuit of softness, SilkySoft


- 3D & Helical Hybrid
- Straight, in-line structure coil design searches for space while entwined with other coils.
- Softest coil in our product history
- Soft coils & the Waving SR makes tight packing & minimal catheter kick-back possible.



Refined Pushability with Unique Gold-colored Non-Slip Coating



Complex Infini SilkySoft

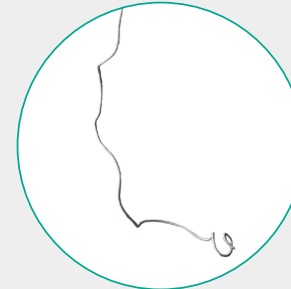

New concept filling coils

- 12 Coil Type SilkySoft
- A soft coil with a 0.012" Primary Coil Diameter for volume maximization.
- Filling coil that moves into space with ease
- Large loops that can fold tightly, configure according to the size of the space and pack them.
- Greater volume gain performance with a single coil
- Due to the Waving SR and softness of the SilkySoft, the coil can fold compactly enabling choice of longer coils possible.



ELECTRO DETACH GENERATOR v4

Sound & Light System To Signal Best Detach Point



EDGv4 indicates by Sound and Light to signal the best Detach Point even when position of 2nd marker of Microcatheter is not detectable under X-Ray Fluoroscopy.